

Technical Data for Alicat **BIO-Series** Mass Flow Meters

10 sccm of Full Scale through 20 slpm of Full Scale



Tel: 888-290-6060

www.alicat.com/bio

Standard Specifications (Contact Alicat for available options.)

SENSOR PERFORMANCE		
Mass Flow Accuracy at calibration conditions ¹	± 0.6% of Reading	16.7% - 100% of Full Scale Range
	± 0.1% of Full Scale	0% - 16.7% of Full Scale Range
Repeatability (2σ)	± (0.1% of Reading + 0.02% of Full Scale)	
Flow Measurement Range	0.01% - 100% of full scale	
Temperature Sensitivity	Mass Flow Zero Shift: ± 0.01% of Full Scale per °C from tare temperature Mass Flow Span Shift: ± 0.01% of Reading per °C from 25°C	
Pressure Sensitivity	Mass Flow Zero Shift: ± 0.01% of Full Scale per atm from tare pressure Mass Flow Span Shift: ± 0.1% of Reading / atm from calibration pressure	
Operating Temperature Range	-10 to 60°C (consult Alicat for expanded range)	
Temperature Accuracy	± 0.75°C	
Operating Pressure Full Scale	160 PSIA (consult Alicat for additional options)	
Pressure Accuracy	Above 1 atm: ± 0.5% of Reading	Below 1 atm: ± 0.07 PSIA
Typical Sensor Response Time	< 10 ms (Adjustable)	
Typical Warm-Up Time	< 1 s	

1 Stated accuracy is after tare under equilibrium conditions. Extreme gas behavior (especially near state boundaries) can introduce additional flow uncertainties. Consult Alicat if higher accuracy is required.

MECHANICAL	
Minimum Operating Pressure	11.5 PSIA common mode pressure (consult Alicat for lower operating pressures) Differential pressure must exceed model pressure drop, see below for details
Maximum Operating Pressure	Damage possible above 175 PSIA common mode pressure Damage possible above 75 PSID differential pressure
Ingress Protection	IP40 (consult Alicat for additional options)
Humidity Range	0 to 95% non-condensing
Dimensions, pressure drop, weight, and process connection specifications are listed on mechanical drawing pages	

CONTROL AND COMMUNICATIONS		
Analog I/O	0-5 VDC Serial and Modbus RTU only	
Digital I/O Options	DeviceNet, EtherCAT, EtherNet/IP, Modbus RTU over RS-232, Modbus RTU over RS-485, Modbus TCP/IP, Profibus, RS-232 Serial, RS-485 Serial	
Electrical Connection	8 pin M12 or Protocol Dependent	
Power Requirements ²	12-24 VDC, 325 mA min.	
Data Update Rate ²	Serial: 40 Hz at 19200 baud	Analog: 1000 Hz
Display Update Rate	10 Hz	
Analog Signal Accuracy	± 0.1% of Full Scale additional uncertainty	

2 Consult the individual operating bulletins for specific industrial protocol power requirements and data transmission specifications.


FEATURES	
STP Reference Conditions	25°C and 1 atm (Default), user configurable
NTP Reference Conditions	0°C and 1 atm (Default), user configurable
Color TFT Display with integrated touchpad	Simultaneously displays Mass Flow, Volumetric Flow, Pressure and Temperature
Gas Select™	98 user selectable gases stored internally. Each gas optimized to match NIST's REFPROP 10 gas property calculations across the operating temperature and pressure ranges for highest accuracy.
COMPOSER™	Allows 20 user definable gas mixes Up to 5 constituent gases per mix, down to percentages of 0.01%

Wetted Materials

FLOW BODY WETTED MATERIALS	OPTION	SENSOR WETTED MATERIALS
316L Stainless Steel, USP VI FDA Certified Viton Elastomers	A	316L Stainless Steel
ASME BPE-2016 Compliance Requires Sensor A Each meter has 2 parts: Flow body and Sensor	B	Polyamide, Alumina, Ceramic, Glass, Gold, Silicon, Nylon, Delrin, Heat Cured Epoxy, RTV, Silicone

BIO-Series

- 0 – 10 sccm
- 0 – 20 sccm
- 0 – 50 sccm

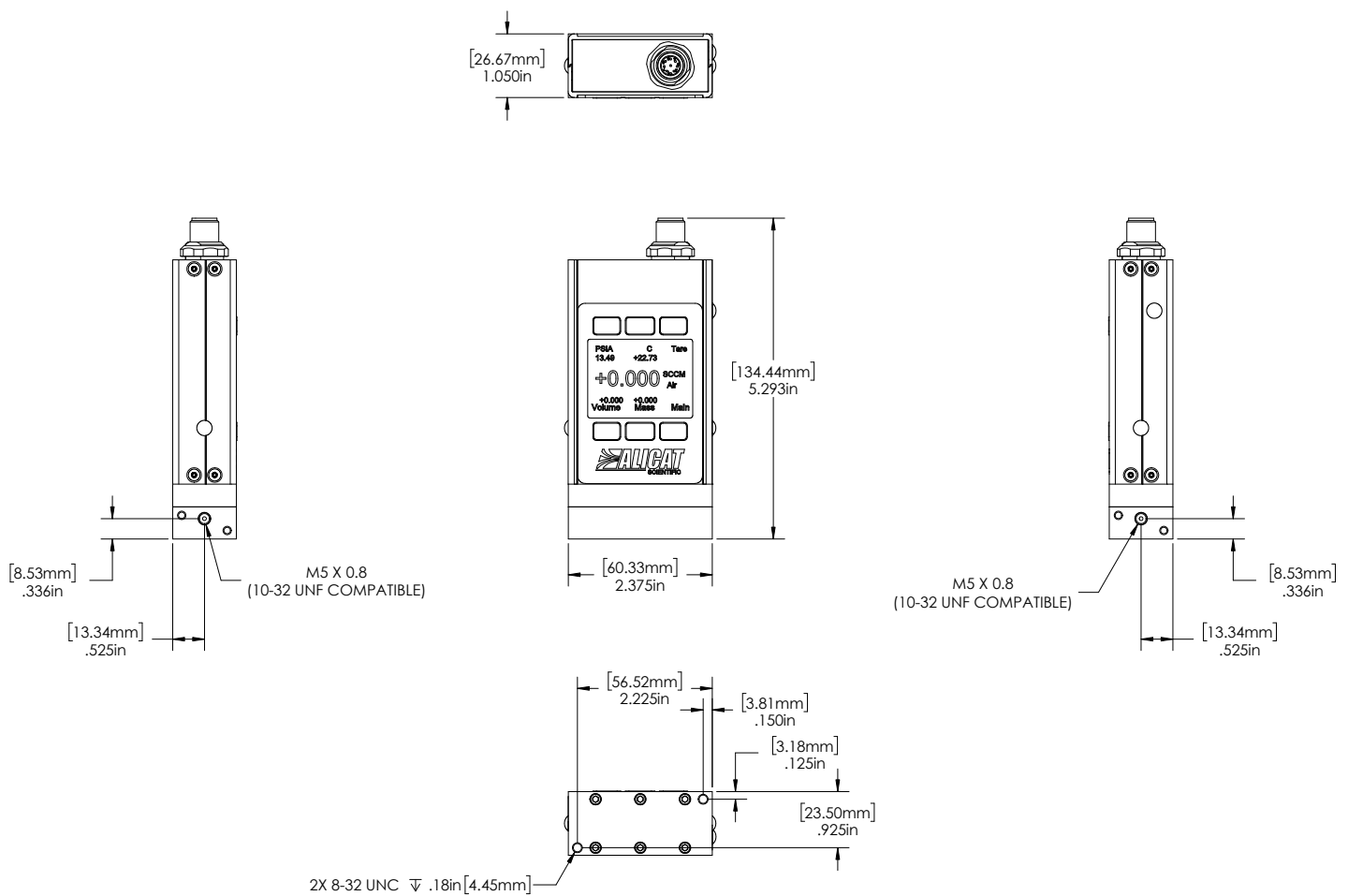

Flow Range Specific Specifications

FULL SCALE FLOW MASS METER	PRESSURE DROP AT FS FLOW (PSID) VENTING TO ATMOSPHERE ³	PHYSICAL DIMENSIONS ⁴	APPROXIMATE WEIGHT	PROCESS CONNECTIONS ⁵
10 sccm to 50 sccm	1.0	5.3"H x 2.4"W x 1.1"D	0.8 lb	M-5 (10-32) Female ⁶


3 Lower pressure drops available, please see our WHISPER-Series mass flow controllers at www.alicat.com/whisper.

4 See drawings for metric equivalents.

5 Additional process connections available on request. Consult Alicat for more information.

BIO-Series

0 – 100 sccm	0 – 1 slpm
0 – 200 sccm	0 – 2 slpm
0 – 500 sccm	0 – 5 slpm
	0 – 10 slpm
	0 – 20 slpm

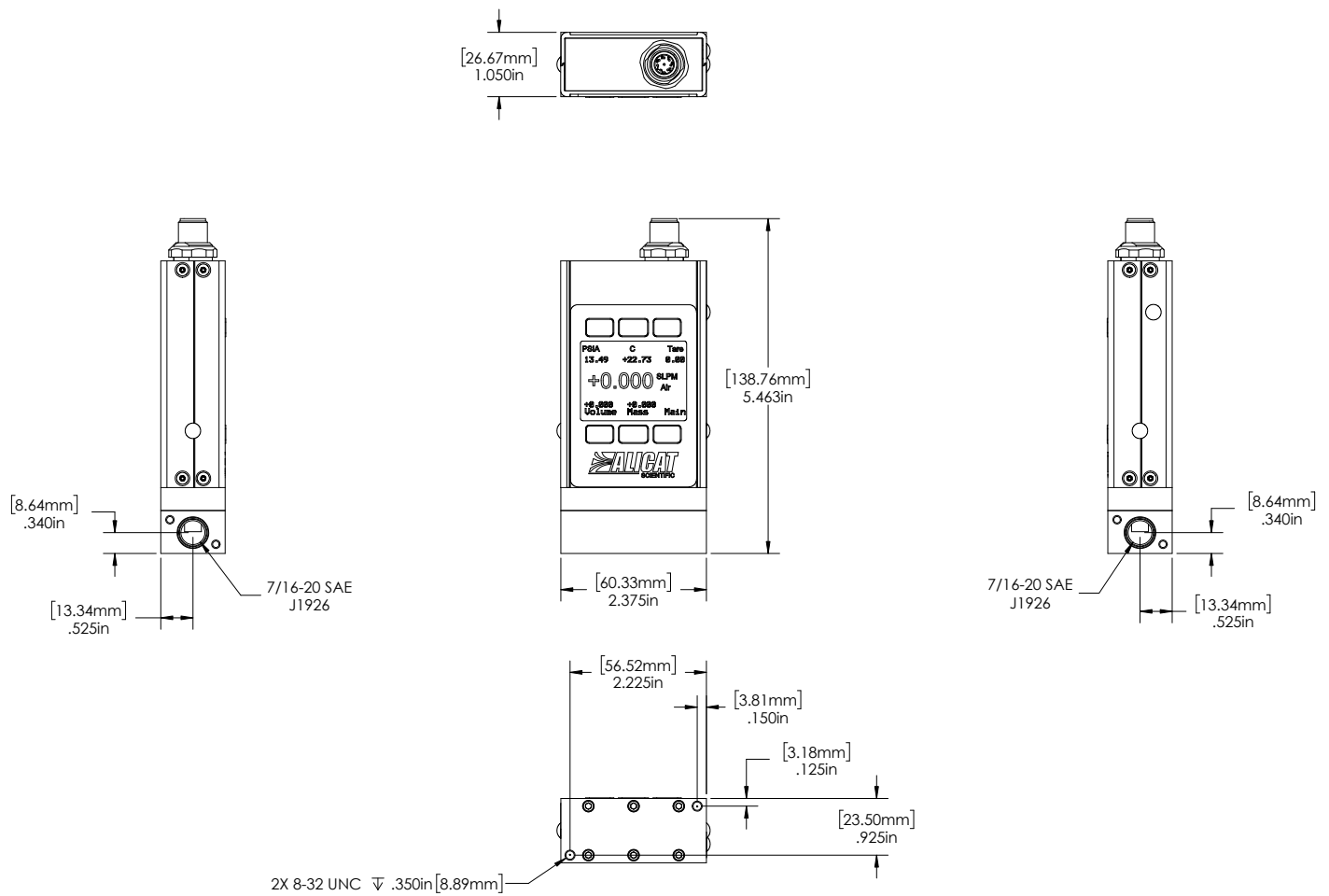

Flow Range Specific Specifications

FULL SCALE FLOW MASS METER	PRESSURE DROP AT FS FLOW (PSID) VENTING TO ATMOSPHERE ³	PHYSICAL DIMENSIONS ⁴	APPROXIMATE WEIGHT	PROCESS CONNECTIONS ⁵
100 sccm to 20 slpm	1.0	5.5"H x 2.4"W x 1.1"D	1.0 lb	7/16-20 SAE ⁴

3 Lower pressure drops available, please see our WHISPER-Series mass flow controllers at www.alicat.com/whisper.

4 See drawings for metric equivalents.

5 Additional process connections available on request. Consult Alicat for more information.

